

 Orthofix Inc.  
3451 Plano Parkway  
Lewisville, Texas 75056-9453 U.S.A.  
1-214-937-3199  
1-888-298-5700  
www.orthofix.com

Australian Sponsor  
Emergo Australia  
Level 20, Tower II  
Darling Park  
201 Sussex Street  
Sydney, NSW 2000  
Australia

 Medical Device Safety Service  
(MDSS)  
Schiffgraben 41  
30175 Hannover  
Germany  
+49 511 6262 8630  
www.mdss.com

---

Device System Name:  
**Firebird® Spinal Fixation System**  
which includes:  
**Firebird® System**  
**Firebird® Deformity System**  
**Firebird® NXG Spinal Fixation System**  
**Phoenix® Minimally Invasive Spinal Fixation System**  
**Phoenix® CDX™ Minimally Invasive Spinal Fixation System**  
**JANUS® Midline Fixation Screw**  
**JANUS® Fenestrated Screw**

English **EN**

2-4

**Orthofix Inc.**  
 3451 Plano Parkway  
 Lewisville, Texas 75056-9453 U.S.A.  
 1-214-937-3199  
 1-888-298-5700  
 www.orthofix.com

**Australian Sponsor**  
 Emergo Australia  
 Level 20, Tower II  
 Darling Park  
 201 Sussex Street  
 Sydney, NSW 2000  
 Australia

**EC REP** Medical Device Safety Service  
 (MDSS)  
 Schiffgraben 41  
 30175 Hannover  
 Germany  
 +49 511 6262 8630  
 www.mdss.com

## English EN

**Device System Name:**

### Firebird® Spinal Fixation System

which includes:

- Firebird® System**
- Firebird® Deformity System**
- Firebird® NXG Spinal Fixation System**
- Phoenix® Minimally Invasive Spinal Fixation System**
- Phoenix® CDX™ Minimally Invasive Spinal Fixation System**
- JANUS® Midline Fixation Screw**
- JANUS® Fenestrated Screw**

**Description:**

The Firebird Spinal Fixation Systems include temporary, multiple component systems comprised of a variety of non-sterile and sterile single use components made of titanium alloy or cobalt chrome alloy that allow the surgeon to build a spinal implant construct. The systems are attached to the vertebral body and ilium by means of screw or hook fixation to the non-cervical spine. The systems consist of an assortment of rods, multi-axial and mono-axial pedicle screws, set screws, lateral offsets, bone screws, screw bodies, hooks, iliac connectors and sterile packed HA coated bone screws.

A subset of the systems' components may be used in pediatric patients. These components consist of a variety of screws ranging in diameters from 4.0mm to 7.5mm and lengths ranging from 25mm to 60mm.

**Indications for Use:**

The Firebird Spinal Fixation Systems are intended to provide immobilization and stabilization of spinal segments in skeletally mature patients as an adjunct to fusion as a pedicle screw fixation system (T1-S2/Ilium), in the treatment of the following acute and chronic instabilities or deformities:

1. Degenerative disc disease (defined as discogenic back pain with degeneration of the disc confirmed by history and radiographic studies).
2. Spondylolisthesis.
3. Trauma (i.e., fracture or dislocation).
4. Spinal stenosis.
5. Deformities or curvatures (i.e., scoliosis, kyphosis, and/or lordosis).
6. Tumor.
7. Pseudoarthrosis.
8. Failed previous fusion.

When used for fixation to the ilium, the offset connectors of the Firebird Spinal Fixation System must be used in conjunction with pedicle screws placed at the S1 or S2 spinal level.

The Firebird Spinal Fixation Systems components are used with certain components of the Spinal Fixation System (SFS), including rods, rod connectors and cross-connectors.

When used for posterior pedicle screw fixation in pediatric patients, the Firebird Spinal Fixation Systems are indicated as an adjunct to fusion to treat adolescent idiopathic scoliosis. Pediatric pedicle screw fixation is limited to a posterior approach.

The Firebird Spinal Fixation Systems are intended to be used with autograft or allograft.

The **Phoenix MIS Fixation System** when used with the Firebird Spinal Fixation Systems is indicated to provide the surgeon with a minimally invasive approach for posterior spinal surgery.

The **JANUS Midline Fixation Screw** and the **JANUS Fenestrated Screw** when used with the Firebird Spinal Fixation Systems is indicated to provide the surgeon with an open, minimally invasive or midline approach for posterior spinal surgery. The JANUS Fenestrated Screws are intended to be used with saline and radiopaque dye.

**Contraindications:**

Contraindications include, but are not limited to:

1. Morbid obesity.
2. Mental Illness.
3. Alcoholism or drug abuse.
4. Pregnancy.
5. Metal sensitivity/allergies.
6. Severe osteopenia.
7. Patients unwilling or unable to follow post-operative care instructions.
8. Use of the Firebird offset connectors for fixation to the ilium is contraindicated when the sacrum is absent or insufficient for implantation of pedicle screws at the S1 or S2 spinal level.
9. Any circumstances not listed under the heading indications.

**Potential Adverse Events:**

All of the possible adverse events associated with spinal fusion surgery without instrumentation are possible. With instrumentation, a listing of possible adverse events includes, but is not limited to:

1. Inability to use pedicle screw fixation due to anatomic limitations (pedicle dimensions, distorted anatomy).
2. Pedicle screw malpositioning, with or without neurological or vascular injury.
3. Proximal or distal junctional kyphosis.
4. Pancreatitis.
5. Pedicle screw failure, such as screw or rod bending, breakage, or loosening, may also occur in pediatric patients, and pediatric patients may be at increased risk for device-related injury because of their smaller stature.
6. Device component fracture.
7. Loss of fixation.
8. Non-union.
9. Fracture of the vertebra.
10. Neurological injury.
11. Vascular or visceral injury.
12. Early or late loosening of any or all of the components.
13. Disassembly and/or bending of any or all components.
14. Foreign body (allergic) reaction to implants, debris, corrosion products, and graft material, including metallosis, straining, tumor formation, and/or auto-immune disease.
15. Pressure on the skin from component parts in patients with inadequate tissue coverage over the implant possibly causing skin penetration, irritation, and/or pain.
16. Post-operative change in spinal curvature, loss of correction, height, and/or reduction.
17. Infection.
18. Pain, discomfort, or abnormal sensations due to the presence of the device.
19. Hemorrhage.
20. Cessation of any potential growth of the operated portion of the spine.
21. Death.

**Note:** Potential risks identified with the use of the device system may require additional surgery.

**Warnings and Precautions:**

1. The safety and effectiveness of this device has not been established for use as part of a growing rod construct. This device is only intended to be used when definitive fusion is being performed at all instrumented levels.
2. The use of pedicle screw fixation in the pediatric population may present additional risks when patients are of smaller stature and skeletally immature. Pediatric patients may have smaller spinal structures (pedicle diameter or length) that may preclude the use of pedicle screws or increase the risk of pedicle screw malpositioning and neurological or vascular injury. Patients who are not skeletally mature undergoing spinal fusion procedures may have reduced longitudinal spinal growth or may be at risk for rotational spinal deformities (the "crankshaft phenomenon") due to continued differential growth of the anterior spine.
3. The implantation of pedicle screw spinal systems in pediatric patients should be performed only by experienced spinal surgeons with specific training in the use of this pedicle screw spinal system in pediatric patients because this is a technically demanding procedure presenting a risk of serious injury to the patient.
4. Preoperative and operating procedures, including knowledge of surgical techniques, good reduction, and proper selection of placement of the implants are important considerations in the successful utilization of the system in pediatric patients.
5. The selection of the proper size, shape and design of the implant for each patient is crucial to the safe use of this device in pediatric patients.
6. The safety and effectiveness of pedicle screw systems have been established only for spinal conditions with significant mechanical instability or deformity requiring fusion with instrumentation. These conditions are: significant mechanical instability or deformity of

the thoracic, lumbar, and sacral spine secondary to severe spondylolisthesis (grades 3 and 4) of the L5-S1 vertebra, degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor, and failed previous fusion (pseudarthrosis). The safety and effectiveness of these devices for any other condition are unknown.

7. The benefit of spinal fusion utilizing any pedicle screw fixation system has not been adequately established in patients with stable spines.
8. Potential risks identified with the use of this device system which may require additional surgery include: device component fracture, loss of fixation, non-union, fracture of the vertebra, neurological injury and vascular or visceral injury.
9. Single use only.
10. Non-sterile; the screws, hooks, rods, dominoes, lateral offsets, spacers, staples, washers, locking nuts, cross connectors, and instruments are sold non-sterile and therefore must be sterilized before use.
11. To facilitate fusion, a sufficient quantity of autologous bone or other appropriate material should be used.
12. Failure to achieve arthrodesis will result in eventual loosening and failure of the device construct.
13. Excessive torque applied to the screws may strip the threads in the bone.
14. DO NOT REUSE IMPLANTS. Discard used, damaged, or otherwise suspect implants.
15. The implantation of pedicle screw spinal systems should be performed only by experienced spinal surgeons with specific training in the use of this pedicle screw spinal system because this is a technically demanding procedure presenting a risk of serious injury to the patient.
16. Based on fatigue testing results, the physician/surgeon should consider the levels of implantation, patient weight, patient activity level, other patient conditions, etc. which may impact the performance of the system.
17. Mixing of dissimilar metals can accelerate the corrosion process. Do not use the titanium alloy or cobalt chrome alloy components of this system with implants of other material composition or components from different manufacturers unless specifically stated.
18. The Firebird Spinal Fixation System and Phoenix MIS Fixation System have not been evaluated for safety and compatibility in the MR environment, nor have the Firebird Spinal Fixation System or the Phoenix MIS Fixation System been tested for heating or migration in the MR environment.
19. Do not attempt to re-sterilize single-use implants that come in contact with body fluids.
20. When using the offset connectors to connect the Firebird spinal construct to the ilium, pedicle screws must be used at the S1 or S2 level of the spine. Do not use the offset connectors to connect the ilium without this intermediate screw fixation.
21. The safety, efficacy and performance of the system have been established for conditions in which the system is used as intended and when used as identified in the Indications for Use. Performance of the system has not been evaluated for use that is contrary to the Intended Use, Indications for Use or for use that is contraindicated. Failure to use the system as indicated could detrimentally affect the performance of its components.
22. Other adverse effects related to pedicle screw fixation, such as screw or rod bending, breakage, or loosening, may also occur in pediatric patients. Pediatric patients may be at increased risk for device-related injury because of their small stature.
23. The correct handling of the implant is extremely important. Implants should not be excessively or repeatedly bent, notched or scratched. These operations can produce defects in surface finish and internal stress concentrations, which may become the focal point for eventual failure of the device.
24. HA coated screws are provided STERILE. Do not use if the package is opened or damaged or if the expiration date has passed.
25. DO NOT re-sterilize the HA coated screws as this could result in injury or require reoperation due to breakage.
26. The safety and effectiveness of this device has not been established when used in conjunction with bone cement or for use in patients with poor bone quality (e.g., osteoporosis, osteopenia). This device is intended only to be used with saline or radiopaque dye.

#### **MRI Compatibility Information:**

The Systems have not been evaluated for safety and compatibility in the MR environment. It has not been tested for heating, migration, or image artifact in the MR environment. The safety of the Systems in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.

#### **Cleaning:**

The HA coated screw implants are sterilized using gamma radiation sterilization. Do not re-sterilize. All other system implants are provided clean but not sterile. Once an implant comes in contact with any human tissue or bodily fluid it should not be re-sterilized or used. Please discard all contaminated implants.

#### **For Firebird Spinal Fixation System Cases 44-9010, 44-9011, 44-9012, 44-9013, 44-9020, 44-9030, 44-9040, 44-9050 and 61-9060:**

All instruments and implants must first be cleaned using established hospital methods before sterilization and introduction into a sterile field. Additionally, all instruments and implants that have been previously taken into a sterile surgical field must first be cleaned using established hospital methods before sterilization and reintroduction into a sterile surgical field. Cleaning can include the use of neutral cleaners followed by a deionized water rinse. All products should be treated with care. Improper use or handling may lead to damage and possible improper functioning of the device.

#### **For All Other Firebird System Cases and Caddies:**

All instruments must be thoroughly cleaned after each use. Cleaning may be done using validated hospital methods or following the validated cleaning processes described below.

#### **From Point of Use:**

Whenever possible, do not allow blood, debris or body fluids to dry on instruments. For best results and to prolong the life of the surgical instrument, reprocess immediately after use.

#### **Instructions for Disassembly and Assembly:**

Prior to cleaning please see the operative technique for disassembly and assembly instructions for the five instruments which require disassembly prior to cleaning, the: Multi-Axial Screw Driver (20-0200), Mono-Axial Screw Driver (20-0300), Multi-Axial Screw Driver (36-1831), Modular Screw Driver (36-1832) and the Midline Modular Screw Driver (36-1833). No other instruments within the system require disassembly prior to cleaning.

#### **Preparation for Cleaning:**

1. Any instruments with moving parts (i.e., knobs, triggers, hinges) should be separated and placed in the open position to allow better access of the cleaning fluid to the difficult to clean areas. Use a soft cloth or plastic bristle brush to remove any visible soil from the instruments.
2. Soak the instruments for a minimum of 10 minutes in sterile water prior to the manual or automated cleaning process.
3. Place scrubbed instruments into the washer baskets.

#### **Manual Cleaning:**

1. After disassembly of the instruments, use a soft cloth and/or a soft plastic bristle brush to remove any visible soil from the outside and inside of the instrument.
2. Prepare Vesphene® Ilse agent at the use-dilution recommended on the label directions (1 ounce per gallon) as follows: Add 1 mL of Vesphene® Ilse to 128.0 mL of potable tap water per the manufacturer's recommendations.
3. Bathe instruments in prepared room temperature solution as recommended by the detergent manufacturer.
4. Manually agitate instruments in Vesphene® solution for 15 minutes.
5. Scrub instruments with a soft plastic bristle brush if visible soil is noted and use a soft plastic bristle brush or pipe cleaner long enough to reach the entire length of the interior lumen to remove the soil.
6. Rinse the instruments in USP <1231> purified water for 1.5 minutes.
7. Hang dry the device.
8. Visually inspect the instruments for visible soil.
9. If visible soil is noted, repeat the steps listed above.

#### **Automated Cleaning:**

1. After disassembly of the instruments use a soft cloth or soft plastic bristle brush to remove any visible soil from the outside of and inside of the instruments. Use a soft plastic bristle brush or pipe cleaner long enough to reach the entire length of the interior lumen to remove the soil from the inner lumens prior to automated cleaning.
2. Set up the washer/disinfectant detergent dose at ½ ounce of Endozime AW Plus® per gallon of water or according to the manufacturer recommendations.
3. Place scrubbed instruments into the washer baskets.
4. Orient instruments in the automated washer's carriers as recommended by the washer manufacturer.
5. The following automated cleaning cycle is recommended (minimum recommended times are provided for each stage):
  - a. Pre-Wash 1: cold potable water, 2 minutes
  - b. Enzyme/Detergent treatment:
    1. Spray, 20 seconds
    2. Soak, 1 minute
    3. Rinse cold potable water, 15 seconds
    4. Rinse cold potable water, 15 seconds
  - c. Wash ≥ 65°C, 2 minutes using Endozime AW Plus®
  - d. Rinse 1: hot potable water, 15 seconds
  - e. Rinse 2: hot potable water, 15 seconds
  - f. Rinse 3: hot potable water, 15 seconds
  - g. Rinse 4: hot potable water, 15 seconds
  - h. Thermal rinse ≥ 93°C for 1 minute
  - i. Heated USP Purified Water Rinse 1: re-circulating 10 seconds
  - j. Heated USP Purified Water Rinse 2: non re-circulating 10 seconds
  - k. Drying: 7 minutes, 115°C
  - l. Visually inspect the instruments for visible soil.
  - m. If visible soil is noted, repeat the above listed steps until no visible soil is noted.

**Note:** Certain cleaning solutions such as those containing caustic soda, formalin, glutaraldehyde, bleach and/or other alkaline cleaners may damage instruments. These solutions should not be used.

**Note:** Visually inspect instruments after cleaning and prior to each use. Discard or return to Orthofix any instruments that are broken, discolored, corroded, have cracked components, pits, gauges, or are otherwise found defective. Do not use defective instruments.

#### **Sterilization:**

The System HA coated screw implants are sterilized using gamma radiation sterilization. Do not re-sterilize. All other implants and instruments are supplied NON-STERILE.

#### **For Firebird Spinal Fixation System Cases 44-9010, 44-9011, 44-9012, 44-9013, 44-9020, 44-9030, 44-9040, 44-9050 and 61-9060:**

The Firebird Spinal Fixation System should be sterilized by the hospital using one of the following recommended cycles when utilizing an FDA cleared sterilization wrap:

Method: Steam	or:	Method: Steam
Cycle: Gravity		Cycle: Prevac
Temperature: 250°F (121°C)		Temperature: 270°F (132°C)
Exposure time: 30 minutes		Exposure time: 8 minutes

#### **For All Other Firebird System Cases and Caddies:**

Prior to use, all implants and instruments should be placed in the appropriate Orthofix case which will be wrapped in an FDA cleared sterilization wrap, or individually wrapped, and placed in the autoclave for sterilization by the hospital using one of the following recommended cycles:

Method: Steam	or:	Method: Steam
Cycle: Gravity		Cycle: Prevac
Temperature: 270°F (132°C)		Preconditioning: Per manufacturer's settings
Exposure time: 15 minutes		Exposure time: 4 minutes
Drying time: 30 minutes		Drying time: 30 minutes
Double wrapped		Double wrapped

#### **Firebird NXG, JANUS Midline Fixation Screw and JANUS Fenestrated Screw Sterilization in Rigid Sterilization Containers:**

When using rigid sterilization containers, clean, inspect and prepare the rigid sterilization container according to the manufacturer's instructions.

Select the appropriate rigid sterilization container with either filtered bottom or solid bottom to properly enclose the Orthofix case(s) (recommended 23¼" long x 11¼" wide container). The following sterilization cycle has been validated:

- Sterilization Method: Steam
- Cycle: Prevac
- Temperature: 270°F (132°C)
- Preconditioning: Per manufacturer's settings
- Exposure time: 4 minutes
- Drying time: 30 minutes

**Patient Information:**

The temporary internal fixation devices used in your recent spinal surgery are metallic implants that attach to the bone and aid in the healing of bone grafts. These implants have been shown to be valuable aids to surgeons in the treatment of bony fusions. These devices do not have the capabilities of living bone. Intact living bone is self-repairing, flexible and occasionally breaks and/or degrades. The anatomy of the human body places a size limitation on any artificial fixation device used in surgery. The maximum size limitation increases the chances of the mechanical complications of loosening, bending or breaking of the devices. Any of these complications could result in the need for additional surgery. Accordingly, it is very important that you follow the recommendations of your physician. Use braces as instructed. By following these instructions, you can increase your chances of a successful result and reduce your risk of injury and/or additional surgery.

**Packaging:**

Packages for each of the components should be intact upon receipt. If a consignment system is used, all sets should be carefully checked for completeness and all components should be carefully checked for damage prior to use. Damaged packages or products should not be used and should be returned to Orthofix.

The System instruments and implants are provided in a modular case specifically intended to contain and organize the system components. The system instruments are organized into trays within the modular case for easy retrieval during surgery. These trays also provide protection to the system components during shipping. Additionally, individual instruments and implants will be provided in sealed poly bags with individual product labels attached to them.

**Note:** The HA coated screws are provided sterile. Do not use if the package is opened or damaged, or if the expiration date has passed.

**Product Complaints:**

Any Health Care Professional (e.g., customer or user of this system of products) who has any complaints or who has experienced any dissatisfaction in the product quality, identity, durability, reliability, safety, effectiveness and/or performance, should notify Orthofix Inc., 3451 Plano Parkway, Lewisville, TX 75056, USA, by telephone at 1-214-937-3199 or 1-888-298-5700 or by e-mail at [complaints@orthofix.com](mailto:complaints@orthofix.com).





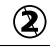







**Further Information:**

Recommended operative techniques for the use of these systems are available upon request from Orthofix at the phone numbers provided above.

**Latex Information:**

The implants, instruments and/or packaging material for the System are not formulated with and do not contain natural rubber. The term "natural rubber" includes natural rubber latex, dry natural rubber, and synthetic latex or synthetic rubber that contains natural rubber in its formulation.

**Caution:** Federal law (U.S.A.) restricts this device to sale by or on the order of a physician.

<b>Rx Only</b>		Federal (U.S.A.) law restricts this device to sale by or on the order of a physician	
	See Instructions for Use		Use By Date
	Orthofix.com/IFU		Manufacturer
	Single Use Only Do Not Reuse		Authorized Representative
	Catalogue Number		Do Not Resterilize
	Provided Non-Sterile		Serial Number
	Sterilized Using Irradiation		Lot Number